

ASSAY FOR THE SECRETED ALKALINE PHOSPHATASE (SEAP)

- 1> - From the 48-hour cotransfected cell cultures, remove 250 μ l of each culture supernatant and transfer into Eppendorf tubes. To be on the safe side, maintain the cultures at 37°C until satisfactory data have been obtained.
- 2> - Heat samples at 65°C for 5 minutes to inactivate endogenous phosphatases (Seap is stable at 65°C). *homoarginine also inactivate endogenous phosphatases*
- 3> - Centrifuge for 2 minutes at room temperature in a Microfuge.
- 4> - Transfer ^{100 μ l} supernatants to new Eppendorf tubes. At this stage, samples may be stored at -20°C indefinitely.
- 5> - In an Eppendorf tube, add 100 μ l of 2 x Seap buffer to 100 μ l of sample. As a Zero standard make up a mixture in triplicate substituting sample with water. Mix on a Vortex.
- 6> - Transfer the contents of each tube to a well of a flat bottom microtiterplate. Avoid creating air bubbles.
- 7> - Incubate plate at 37°C for 10 - 15 minutes.
- 8> - During this incubation make up the p-nitrophenylphosphate solution (Seap enzyme substrate) and prewarm it at 37°C for 5 minutes.
- 9> - Add 20 μ l of the substrate solution to each well, preferably using a multipipetter.
- 10> - Using an ELISA microplate reader with an automatic shaker and incubator unit, measure the OD at (405) nm at regular intervals (e.g. every 5 minutes) over 60 minutes while the plate is being incubated at 37°C. (Program a 5-second shaking before any reading).
- 11> - Calculate the levels of Seap activity at a point on the curve when the changes in OD are linear with respect to time (e.g. at 15 - 30 minutes).

***** * *****

Buffer and Chemicals

✓ 2 x Seap Buffer (for 50 ml)

Amount	Stock	Final Conc.
10.51 g*	diethanolamine (100% sol.)	2 M
50 μ l	1 M MgCl ₂	1 mM
226 mg	L-homoarginine	20 mM

(*) Weigh exactly 10.51 g of diethanolamine in a tared beacker. Add distilled water up to 45 ml. Stir to homogenize. Add 50 μ l of 1 M MgCl₂ while stirring. In a separate 15-ml conical tube, dissolve 226 mg of L-homoarginine in 2-3 ml of distilled water. Add the solution to diethanolamine/MgCl₂ solution under constant stirring. Bring up to 50 ml.

120 mM p-nitrophenylphosphate is made in 1 x Seap buffer (fresh).

Amount G (mg) = (120 mM x 371.12 x vol)/1000

Where: vol = [(# wells x 20 μ l) + 100 μ l extra] / 1000

* Make the solution in 1 x Seap buffer (make fresh)

For 51 wells of the microplate (48 samples + 3 blanks), dissolve 50 mg p-nitrophenylphosphate in 1.120 ml of 1 x Seap buffer.

→ 600 μ l stock + 600 μ l H₂O

Chemicals

Name	Cat No	Storage
Diethanolamine (Fluka)	31589	Room T
L-homoarginine hydrochloride (Sigma)	H-1007	4°C
p-nitrophenylphosphate* (Fluka)	71768	4°C

(*) Also known as 4-nitrophenylphosphate disodium salt hexahydrate

0.21805
0.226
" "

SEAP ASSAY SHEET

I - Assay Title:

* Assay #: 72

Date #: [REDACTED]

Investigator's Name: Wen

*Test Compounds: Mal X & Mal X 3Na⁺

transfected at [REDACTED] 11:30

*Concentrations:

II - DNA Transfection: Ratio (2:1) HIV/SEAP:pcTAT

*HIV/Seap: 0.4 µg/well x 40 = 16.0 µg ==> From 0.376% stock: 42.55 µl

*pcTAT: 0.2 µg/well x 40 = 8.0 µg ==> From 0.639% stock: 12.52 µl

*Total DNA (µg) = 24 µg

* Vol. Cellfectin = Total DNA x 6 = 144 µl

* Vol. 150 mM NaCl = Cellfectin (µl) / 0.6 = 240 µl

HIV/SEAP
42.55 µl

pcTAT
12.52 µl

(A)

SALT

(B)

Cellfectin, 144 µl

* Transfection cocktail (µl)/ well:

$$\frac{240 + 240 + 144 + 42.55 + 12.52}{40} = 16.98 \mu l$$

III. Linbro® 24 flat bottom well of 17 mm

A

10			40	
20			60	

B

10			40	
20			60	

Preparation of Drug conc.

(NOTE)

C

C		C Salt
CT		
C DMSO		
CT DMSO		CT Salt

40: 78
60: 97

40: 23
60: 46

REPEATED-READS #: 0025

07:11:06.56 PM

TABLE OF ABSORBANCE VALUES

FILE: FILE 3 TITLE: FILE 3 REPORT

1 2 3 4 5 6 7 8 9 10 11 12

Mal.4 DMSO											
A	0#000-0#002	0#004-0.217-0.215-0.216	0.064	0.080	0.089	0.126	0.133	0.146			
			0.078		40 $\mu\text{g/mL}$	0.135					
B	-0.210-0.214-0.205-0.215-0.215-0.213	0.151	0.153	0.207	0.124	0.153	0.155				
			0.170		0.144						
C	0.070	0.016	0.060	0.044	0.214-0.214	0.105	0.089	0.131	0.362	0.171	0.169
0.641			0.048		20 $\mu\text{g/mL}$	0.168		60 $\mu\text{g/mL}$	0.170		
D	0.689	0.649	0.663	0.756	0.208-0.210	0.126	0.122	0.120	0.362	0.172	0.197
			0.689			0.123		0.185			
C DMSO											
E	0.094	0.100	0.081	0.110	0.212-0.209	0.093	0.067	0.087	0.101	0.080	0.100
0.575			0.096		10 $\mu\text{g/mL}$	0.082		40 $\mu\text{g/mL}$	0.094		
F	0.680	0.711	0.680	0.614	0.211-0.213	0.830	0.776	0.832	0.649	0.587	0.501
			0.671			0.813		0.579			
C NaOH											
G	0.064	0.067	0.080	0.079	0.203-0.211	0.089	0.114	0.100	0.114	0.118	0.098
0.632			0.073		20 $\mu\text{g/mL}$	0.101		60 $\mu\text{g/mL}$	0.110		
H	0.784	0.708	0.643	0.684	0.209-0.205	0.681	0.720	0.609	0.481	0.437	0.437
			0.705			0.670		0.452			

Mal.4 - DMSO Ref: DMSO
% inhibition

10 $\mu\text{g/mL}$: 84
20 $\mu\text{g/mL}$: 97
40 $\mu\text{g/mL}$: 98
60 $\mu\text{g/mL}$: 97

Mal.4 - NaOH Ref: NaOH
% inhibition

10 $\mu\text{g/mL}$: 0 (-16)
20 $\mu\text{g/mL}$: 10
40 $\mu\text{g/mL}$: 23
60 $\mu\text{g/mL}$: 46

SEAP assay of NDGA derivatives

1. Detach CoS7 cells ^{from} ~~the~~ four 90mm culture plates by ^{adding} 0.05% Trypsin (0.05% Trypsin in CMF-PBS, 1mM EDTA) 0.5 ml after 2 times CMF-PBS Wash.
2. Inactivate trypsin by adding 5ml complete ^{IMDM} medium (10% Fetal Bovine Serum, antibiotics) to each plate.
3. Suspend the cells by gentle pipetting. Count by Hemocytometer 160×10^4 cells/cm³ is the cell density.
4. Seed ~~the~~ 80ul cell suspension into 10 Linbro 24-wells culture plates which ^{were} precoated with 0.1% gelatin & contained 0.3ml complete medium each well. (1.3×10^5 cells/Well) (18:00)
5. Transfect the cells after 26h incubation (37°C, 95% Air-5% CO₂) by Adding 30 ul ppt (0.6 μ g DNA/Well) to each well that containing 0.3ml fresh complete medium. The control wells accept same ppt without DNA. (20:00)
Ca-PO₄ ppt preparation:
 - a. for control wells: Ca²⁺ soln (0.25M CaCl₂, 0.025M HEPES, pH 7.1) 1ml drop into 1ml bubbling PO₄³⁻ soln (0.28M NaCl, 0.025M HEPES, ^{0.0015} ~~0.005~~ M Na₂PO₄, pH 7.1).
 - b. for transfect wells:
Three 2-ml tubes containing ~~pBC12/HIV/SEAP 18ug, pBC12/CMV/tz 9ug~~ 675ul PO₄³⁻ soln were added Ca²⁺ soln 675ul ~~con~~ with DNAs (pBC12/HIV/SEAP 18ug, pBC12/CMV/tz 9ug) dropwise with bubbling. then sit 30min before use.
6. Incubate at 37°C for 18h.

7
6 Map of Plates:

I

	C			C #1	20μM
	C Mal.4	20μM		C #2	20μM
	CT			CT #1	20μM
	CT Mal.4	20μM		CT #2	20μM

A # Mal.4

C	0μM	C	10μM	C	60μM
C	Mal.4 3μM	C	30μM	C	100μM
CT	0μM	CT	10μM	CT	60μM
CT	3μM	CT	30μM	CT	100μM

C #2

C	0μM	C	10μM	C	60μM
C	3μM	C	30μM	C	100μM
CT	0μM	CT	10μM	CT	60μM
CT	3μM	CT	30μM	CT	100μM

E #4

C	0μM	C	10μM	C	60μM
C	3μM	C	30μM	C	100μM
CT	0μM	CT	10μM	CT	60μM
CT	3μM	CT	30μM	CT	100μM

G #6

C	0μM	C	10μM	C	60μM
C	30μM	C	30μM	C	100μM
CT	0μM	CT	10μM	CT	60μM
CT	30μM	CT	30μM	CT	100μM

II

	C #3	20μM	C #4	20μM	C #5	20μM	C #6	20μM
	C #7	20μM	CT #3	20μM	CT #4	20μM	CT #5	20μM
	CT #7	20μM						

B #1

C	0μM	C	10μM	C	60μM
C	3μM	C	30μM	C	100μM
CT	0μM	CT	10μM	CT	60μM
CT	3μM	CT	30μM	CT	100μM

D #3

C	0μM	C	10μM	C	60μM
C	3μM	C	30μM	C	100μM
CT	0μM	CT	10μM	CT	60μM
CT	3μM	CT	30μM	CT	100μM

F #5

C	0μM	C	10μM	C	60μM
C	3μM	C	30μM	C	100μM
CT	0μM	CT	10μM	CT	60μM
CT	3μM	CT	30μM	CT	100μM

H #7

C	0μM	C	10μM	C	60μM
C	3μM	C	30μM	C	100μM
CT	0μM	CT	10μM	CT	60μM
CT	3μM	CT	30μM	CT	100μM

C = Control wells (without DNA)
CT = Transfect wells (with DNA)

- 8 Remove growth medium, add 0.500ul medium containing test compounds. The final concentration of DMSO is 0.3% ([REDACTED])

Δ Dilution of test compounds:

Stock: Mono Me: 10.53 $\mu\text{g}/\text{ul}$ USE 3ul in 1ml medium = 100 μM

DiMe: 11.0 $\mu\text{g}/\text{ul}$ "

TriMe: 11.41 $\mu\text{g}/\text{ul}$ "

tetraMe: 11.93 $\mu\text{g}/\text{ul}$ "

For lower concentration, dilute the stock by DMSO, & use 3ul in 1ml medium to keep the final DMSO concentration is 0.3%.

Δ DMSO stock was ~~add~~^{mixed} into medium just prior to use by vortexing.

- 9 Incubate for 48h. Remove 200ul medium from each well. for SEAP assay. ([REDACTED] 14:00).

- 10 The medium used in SEAP assay is 10ul for each sample.

REPEATED READS #1 0025

05:50:03.26 PM

TABLE OF ABSORBANCE VALUES

FILE: FILE 3 TITLE: FILE 3 REPORT

1 2 3 4 5 6 7 8 9 10 11 12

I						II						
A	-0#024	0#005	0#020	-0.217	-0.216	-0.217	0.011	0.012	0.013	0.004	0.005	0.002
							0.012			0.004	0.007	0.001
E	0.012	-0.008	0.012	0.012	0.020	0.008	0.002	0.004	0.005	0.003	0.019	0.003
		0.005			0.013			0.004				
C	-0.004	0.004	0.016	0.004	0.008	-0.001	0.585	0.704	0.631	0.558	0.394	0.421
		0.005			0.004			0.646	#3	0.530	0.411	0.416
D	1.126	1.040	1.102	0.801	0.766	0.811	0.347	0.373	0.326	0.502	0.428	0.410
Mal. %		1.089			0.793	#1		0.349	#7	#4	#5	#6
E	1.010	0.855	0.915	0.890	0.753	0.799	0.018	0.020	0.007	0.012	0.006	0.023
		0.927			0.814	#2	0.019		0.010		0.009	
F	-0.213	-0.216	-0.215	-0.213	-0.213	-0.216	0.008	0.014	0.007	-0.010	0.010	0.036
							0.011		-0.002		0.023	
G	0.007	0.011	0.013	0.005	0.001	0.017	1.323	1.433	1.105	1.111	0.292	0.263
	0.009		0.009		0.009		1.378		1.108		0.281	
H	0.007	-0.001	0.016	0.003	0.009	0.026	1.284	1.234	0.783	0.715	0.165	0.162
	0.003		0.010		0.018		1.559		0.749		0.164	

20 μM
inhibition.

Mal. % 14.9 Mal. %

#1 : 28.0 =

#2 : 25.3 =

#3 : 42.1 =

#4 : 51.5 =

#5 : 62.8 =

#6 : 61.7 =

#7 : 68.2 =

Mal. % % inhibition

0.3 μM : 8.2

10 μM : 19.2

30 μM : 44.7

60 μM : 80.0

100 μM : 89.6

Exhibit B(3)

Page 4

07:46:07.64 PM

TABLE OF ABSORBANCE VALUES

FILE: FILE 3 TITLE: FILE 3 REPORT

1 2 3 4 5 6 7 8 9 10 11 12

D											
A	0#002-0#004	0#003-0.211	-0.209	-0.212	0.008	0.016	0.003	-0.001	0.002	0.016	
					#2	0.012		0.001		0.009	
#2	0.002	0.002	0.012	0.012	0.015	0.006	0.003	0.009	0.000	-0.005	0.010
	0.002		0.012		0.011		0.006		-0.003		0.010
C	-0.005	0.007	0.008	0.004	0.007	-0.001	1.200	1.149	0.871	1.088	0.374
	0.001		0.006		0.003		1.175		0.980		0.366
D	1.385	1.306	1.025	1.098	0.423	0.464	1.067	1.185	0.570	0.530	0.236
	1.346		1.062		0.444		1.126		0.550		0.231
E											
E	1.288	1.327	0.905	1.005	0.330	0.314	0.017	0.012	0.002	0.007	0.003
	1.308		0.955		0.322		#4	0.015		-0.005	0.011
F	-0.202	-0.210	-0.210	-0.210	-0.209	-0.209	0.006	0.010	0.009	0.002	0.004
							0.008		0.006		0.010
B											
G	1.239	1.354	1.209	1.239	0.420	0.350	1.122	1.205	0.921	0.933	0.208
#1	1.297		1.224		0.385		1.164		0.927		0.243
H	1.387	1.279	0.682	0.580	0.203	0.193	1.178	1.165	0.475	0.435	0.085
	1.333		0.631		0.198		1.172		0.453		0.100

#1	#2	#3	#4
% inhibition	% inhibition	% inhibition	% inhibition
3 μM : -3.3(0)	3 μM : 2.8	3 μM : 3.7	3 μM : -1.3(0)
10 " : 5.7	10 " : 11.9	10 μM : 15.8	10 μM : 19.8
30 " : 52.8	30 " : 29.4	30 μM : 52.5	30 μM : 60.9
60 " : 70.8	60 " : 67.8	60 μM : 69.3	60 μM : 82.4
100 " : 86.0	100 " : 76.3	100 μM : 81.0	100 μM : 92.2

REPEATED-READS #: 0025

09:06:41.21 PM

TABLE OF ABSORBANCE VALUES

FILE: FILE 3 TITLE: FILE 3 REPORT

1 2 3 4 5 6 7 8 9 10 11 12

A 0#000-0#008 0#010-0.210-0.208-0.210-0.209-0.209-0.210-0.208-0.209-0.209

B -0.207-0.208-0.209-0.208-0.208-0.209-0.208-0.209-0.209-0.208-0.209-0.211

G #6

H #7

C	0.012	0.014	0.012	0.022	0.017	0.020	0.011	0.014	0.011	0.014	0.019	0.015
	0.013		0.017		0.019		0.013		0.013		0.017	
D	0.007	0.013	0.015	0.014	0.026	0.012	0.010	0.019	0.004	0.020	0.038	0.043
	0.010		0.015		0.019		0.015		0.012		0.041	
E	1.192	1.134	0.699	0.687	0.277	0.262	1.391	1.489	0.658	0.623	0.317	0.300
	1.163		0.693		0.270		1.440		0.641		0.309	
F	1.215	1.138	0.374	0.400	0.202	0.184	1.163	0.974	0.380	0.317	0.189	0.239
	1.177		0.387		0.193		1.069		0.349		0.214	

G -0.205-0.208-0.208-0.206-0.199-0.208-0.206-0.207-0.207-0.207-0.209-0.210

H -0.207-0.207-0.204-0.206-0.206-0.204-0.189-0.207-0.207-0.206-0.206-0.210

#6

% inhibition

#7

% inhibition

3 μM : -1.5
10 μM : 41.0
30 μM : 67.7
60 μM : 78.2
100 μM : 84.9

3 μM : 26.1
10 μM : 56.0
30 μM : 76.4
60 μM : 79.5
100 μM : 87.9

REPEATED-READS #: 0025

10:33:14.31 PM

TABLE OF ABSORBANCE VALUES

FILE: FILE 3 TITLE: FILE 3 REPORT

1 2 3 4 5 6 7 8 9 10 11 12

A -0#005-0#012 0#018-0.209-0.210-0.211-0.210-0.209-0.208-0.209-0.210-0.211

E -0.208-0.209-0.209-0.209-0.209-0.209-0.209-0.209-0.209-0.209-0.207-0.211

#5 B

0.028 0.031 0.037 0.037 0.040 0.041 0.206-0.207-0.209-0.208-0.208-0.209
0.030 0.037 0.041

D 0.015 0.028 0.036 0.028 0.047 0.035 0.208-0.210-0.208-0.208-0.208-0.209
0.022 0.032 0.041

E 1.500 1.394 0.701 0.680 0.318 0.291 0.207-0.207-0.208-0.207-0.209-0.209
1.44 0.691 0.305

F 1.131 1.033 0.385 0.417 0.222 0.185 0.207-0.208-0.207-0.206-0.208-0.212
1.082 0.401 0.504

G -0.206-0.206-0.208-0.207-0.206-0.208-0.206-0.205-0.207-0.207-0.208-0.208

H -0.207-0.204-0.209-0.208-0.207-0.207-0.207-0.207-0.207-0.207-0.207-0.207

#5 % inhibition

3 μM : 25.2

10 μM : ~~53.8~~ 53.8

30 μM : 74.0

60 μM : 81.4

100 μM : 88.5

% Inhibition of NDGA Derivatives on SEAP assay.

	Max. 4	#1	#2	#3	#4	#5	#6	#7
3 μ M	8.2	0 ⁻³³	2.8	3.7	0 ⁻¹³	25.2	0 ⁻¹⁵	26.1
10 μ M	19.2	5.7	21.9	15.8	19.8	53.8	40.1	56.0
30 μ M	44.7	52.8 52.8	29.4	52.5	60.9	74.0	67.7	76.4
60 μ M	80.0	70.8 70.8	67.8	69.3	82.4	81.4	78.2	79.5
100 μ M	89.6	86.0	76.3	81.0	92.2	88.5	84.9	87.9

Figure B(1)

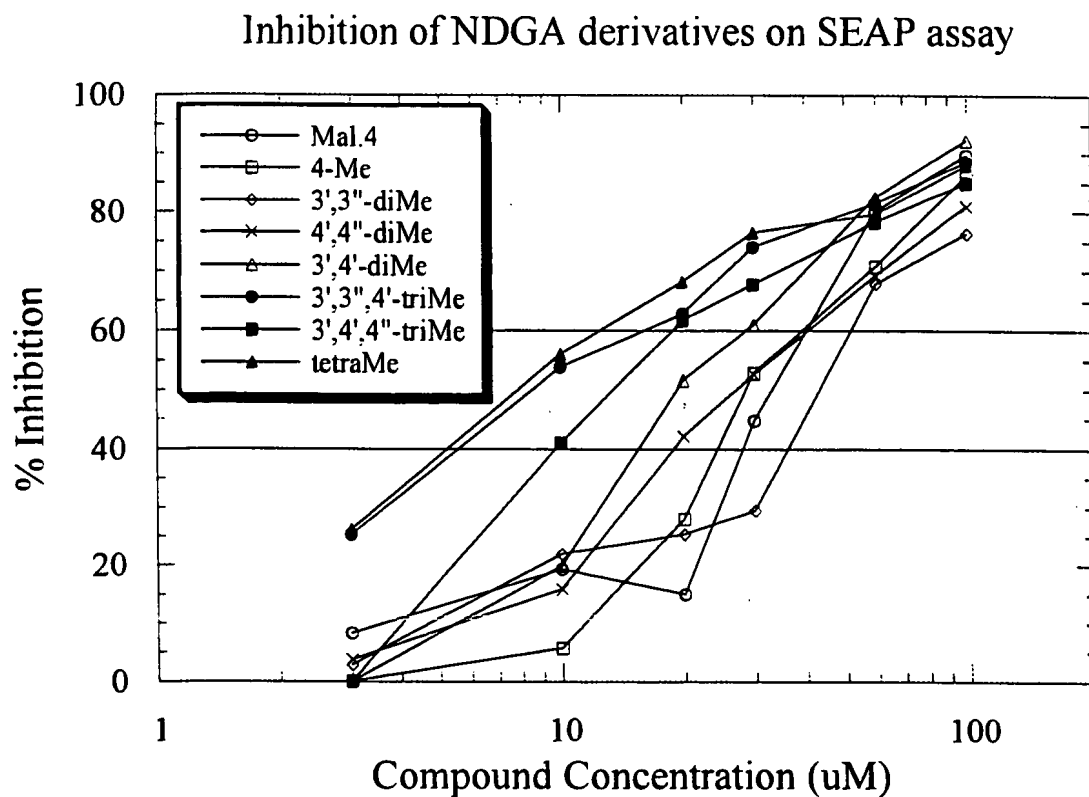


Figure B(2)

